**HOMEDICS**

THE LEADERS IN HEALTH & WELLBEING

# DUO CLIMATE™

HOT AND COLD SONIC FACIAL WAND

PLEASE TAKE A  
MOMENT NOW

REGISTER YOUR PRODUCT AT:

[www.homedics.com/register](http://www.homedics.com/register)

YOUR VALUABLE INPUT REGARDING  
THIS PRODUCT WILL HELP US CREATE  
THE PRODUCTS YOU WILL WANT  
IN THE FUTURE.

INSTRUCTION MANUAL AND WARRANTY INFORMATION

ONE YEAR LIMITED WARRANTY

FAC-420

## **IMPORTANT SAFETY INSTRUCTIONS:**

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:  
**READ ALL INSTRUCTIONS BEFORE USE.**

### **DANGER – TO REDUCE THE RISK OF ELECTRIC SHOCK:**

- DO NOT place or store appliance where it can fall or be pulled into a tub or sink.
- DO NOT place in or drop into water or other liquid.
- DO NOT reach for a device that has fallen into water. Unplug immediately.

### **WARNING – TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK, OR INJURY TO PERSONS:**


- Close supervision is necessary when this appliance is used by, on, or near children, invalids or disabled persons.
- Use this product only for its intended use as described in this manual.
- Keep the cord away from heated surfaces.
- Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
- Never use while sleeping or drowsy.

### **SAVE THESE INSTRUCTIONS**

- This product has no serviceable parts. This product is for household use only. Do not attempt to repair this product.
- If you are under a doctor's care for any skin condition or if you think you may have a skin condition, please consult your physician before using this product.
- Do not use on broken skin or areas of the skin with rash, severe acne, or other skin conditions.
- Individuals with pacemakers and pregnant women should consult a doctor before using this appliance.

This product has an internal, non-replaceable lithium battery. Please dispose in accordance with your local, state, province, and country guidelines.

# DUO CLIMATE HOT AND COLD SONIC FACIAL WAND
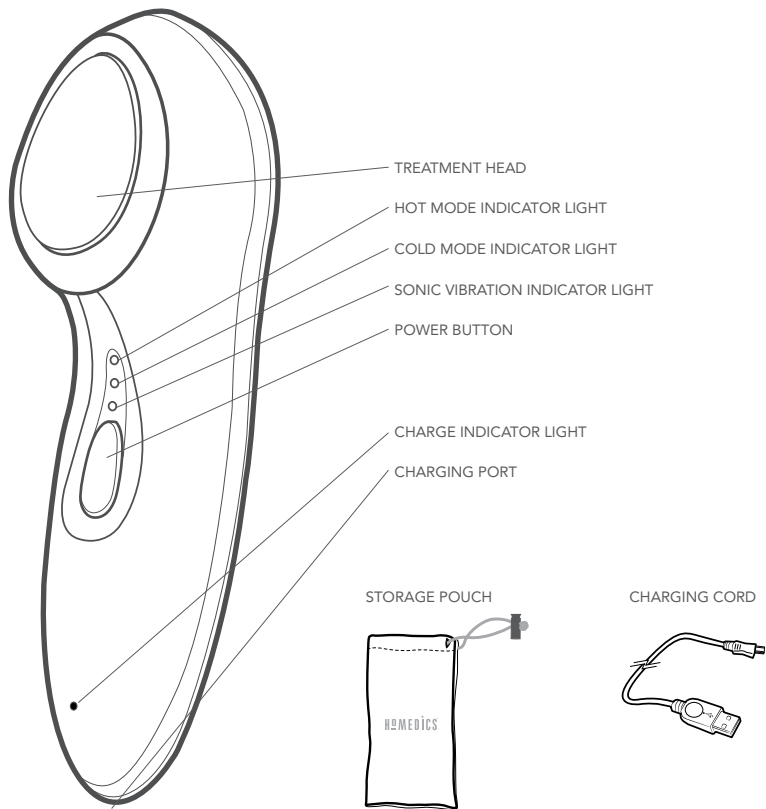
# UNIQUE FEATURES & SPECIFICATIONS

Ever hear the trick of putting eye serums in the fridge? Or rubbing lotion between your palms to warm it up? Temperature dramatically affects how your skin interacts with beauty products. With the DUO Climate Hot and Cold Sonic Facial Wand, gentle heat, refreshing coolness, and sonic massage are just one touch away. Use the hot mode to infuse creams and serums, then switch to the cold mode to tighten and seal skin while reducing redness and puffiness (it's great for under-eye bags!). Plus sonic vibration, which boosts absorption while giving skin a more youthful glow.

## DUAL HOT AND COLD MODES

Supercharges daily cleansing routine

## 4 LEVELS OF TREATMENT

Adjusts for the perfect combination of comfort and effectiveness

## SONIC VIBRATION

Boost effectiveness of Hot and Cold treatment modes

## HOT MODE

Helps to infuse your favorite creams and serums

## COLD MODE

Helps to reduce redness and puffiness

## CORDLESS AND RECHARGEABLE

Charging Cord and Storage Pouch included


## AUTO-TIMER

The auto-timer is programmed to run for 2.5 minutes per setting

# HOW TO USE


## CHARGING THE BATTERY

Before turning the unit on, plug the charging cable into the charging port on the unit and connect to any computer to charge. You can also plug it into a household electrical outlet with a USB adapter plug (not included). The charge indicator light will flash red while it's charging, then turn solid blue when it's fully charged. A full charge takes about 3 hours and lasts for up to 60 minutes of use.
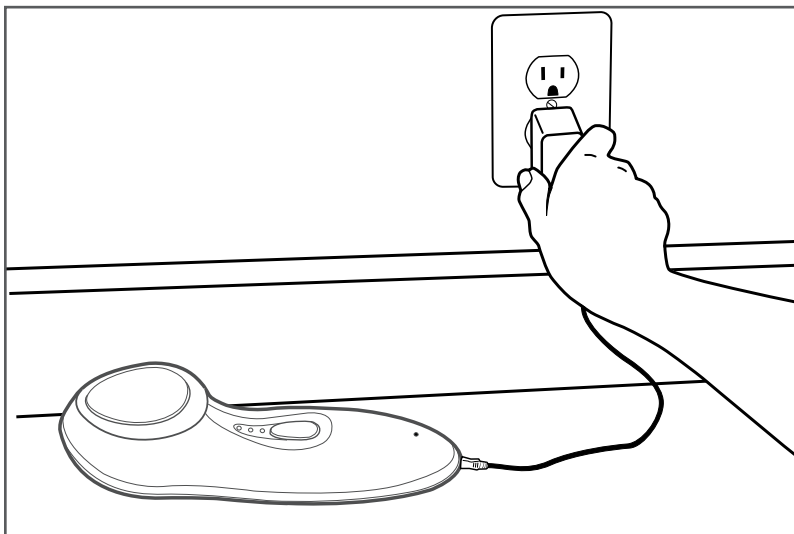
# HOW TO USE


After washing your face, apply your favorite moisturizer, cream or serum.
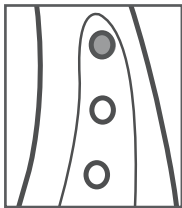
## HOT MODE

- Press the POWER button on DUO Climate Hot and Cold Sonic Facial Wand once to activate HOT MODE. The red indicator light will illuminate. Press the POWER button once more to incorporate soothing Sonic Vibration with HOT MODE. When the Sonic Vibration is activated with HOT MODE, both the red and yellow indicator lights will be illuminated. Glide the wand across your face, using a gentle, circular motion over your nose, cheeks, forehead and chin until your skincare product is absorbed.

**CAUTION:** Do not use hot mode directly on eyelids. Do not hold the DUO Climate Hot and Cold Sonic Facial Wand on the same area of the skin for more than 10 seconds while using heat. The product is designed to move continuously across the skin.
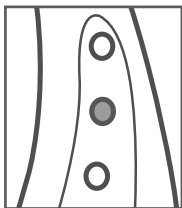
## COLD MODE


- Press the POWER button a third time to activate COLD MODE, which tightens and seals the skin while reducing the appearance of redness and puffiness. The blue indicator light will now only be illuminated. Press the POWER button one more time to incorporate the sonic vibration with COLD MODE. Both the blue and yellow indicator lights will be illuminated. Place the wand on any areas you want to soothe, focusing on problem areas like puffiness under your eyes.

- Press and hold the POWER button to turn the Sonic Facial Wand off.

**NOTE:** For best results, use the Sonic Facial Wand every day after washing your face.

# CLEANING AND CARE

**CAUTION:** before cleaning the unit, turn power off and detach unit from charging cord.
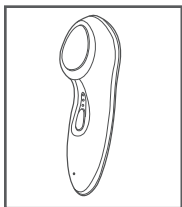
## TO CLEAN THE SURFACE

Clean the surface of the unit with a clean, damp, soft cloth.

**NOTE:** Never use liquids or abrasive cleaners to clean the housing of the unit.

## TO STORE

Place components into the storage pouch and store in a safe, dry place.

# TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
<b>Unit will not operate</b>	<ul style="list-style-type: none"><li>• Unit is not charged</li></ul>	<ul style="list-style-type: none"><li>• Charge the unit, plugging the charging cord into the unit until secure.</li><li>• Press and hold the POWER BUTTON to turn the power on.</li><li>• Check circuits, fuses, try a different outlet.</li></ul>
<b>Unit seems to be losing power</b>	<ul style="list-style-type: none"><li>• Unit is not fully charged</li></ul>	<ul style="list-style-type: none"><li>• Charge the unit, plugging the charging cord into the unit until secure.</li></ul>
<b>Treatment Head is hot when turning the power off after using the cold mode</b>	<ul style="list-style-type: none"><li>• This is normal</li></ul>	<ul style="list-style-type: none"><li>• Allow unit to cool</li></ul>

# FCC STATEMENT

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment. Such modifications could void the user authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This class B digital device complies with Canadian ICES-003.

THE LEADERS IN HEALTH & WELLBEING

## ONE YEAR LIMITED WARRANTY

HoMedics sells its products with the intent that they are free of defects in manufacture and workmanship for a period of one year from the date of original purchase, except as noted below. HoMedics warrants that its products will be free of defects in material and workmanship under normal use and service. This warranty extends only to consumers and does not extend to Retailers.

To obtain warranty service on your HoMedics product, contact a consumer relations representative for assistance. Please make sure to have the model number of the product available.

HoMedics does not authorize anyone, including, but not limited to, Retailers, the subsequent consumer purchaser of the product from a Retailer or remote purchasers, to obligate HoMedics in any way beyond the terms set forth herein. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; improper installation; unauthorized repairs or modifications; improper use of electrical/power supply; loss of power; dropped product; malfunction or damage of an operating part from failure to provide manufacturer's recommended maintenance; transportation damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of HoMedics.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adoption to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. THERE SHALL BE NO OTHER WARRANTIES EXPRESS OR IMPLIED INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS OR ANY OTHER OBLIGATION ON THE PART OF THE COMPANY WITH RESPECT TO PRODUCTS COVERED BY THIS WARRANTY. HOMEDICS SHALL HAVE NO LIABILITY FOR ANY INCIDENTAL, CONSEQUENTIAL OR SPECIAL DAMAGES. IN NO EVENT SHALL THIS WARRANTY REQUIRE MORE THAN THE REPAIR OR REPLACEMENT OF ANY PART OR PARTS WHICH ARE FOUND TO BE DEFECTIVE WITHIN THE EFFECTIVE PERIOD OF THE WARRANTY. NO REFUNDS WILL BE GIVEN. IF REPLACEMENT PARTS FOR DEFECTIVE MATERIALS ARE NOT AVAILABLE, HOMEDICS RESERVES THE RIGHT TO MAKE PRODUCT SUBSTITUTIONS IN LIEU OF REPAIR OR REPLACEMENT.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resaled products, including but not limited to sale of such products on Internet auction sites and/or sales of such products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate as to any products or parts thereof which are repaired, replaced, altered, or modified, without the prior express and written consent of HoMedics.

This warranty provides you with specific legal rights. You may have additional rights which may vary from country to country. Because of individual country regulations, some of the above limitations and exclusions may not apply to you. For more information regarding our product line in the USA, please visit: [www.homedics.com](http://www.homedics.com). For Canada, please visit: [www.homedics.ca](http://www.homedics.ca).

### FOR SERVICE IN THE USA:

EMAIL: [CSERVICE@HOMEDICS.COM](mailto:CSERVICE@HOMEDICS.COM)  
MONDAY-FRIDAY 8:30AM-7:00PM EST  
1.800.466.3342

### FOR SERVICE IN CANADA:

EMAIL: [CSERVICE@HOMEDICSGROUP.CA](mailto:CSERVICE@HOMEDICSGROUP.CA)  
MONDAY-FRIDAY 8:30AM-5:00PM EST  
1.888.225.7378

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©2017 HoMedics, LLC. Tous droits réservés. HoMedics est une marque déposée de HoMedics, LLC.