Well at  
*Walgreens*

PREMIUM WRIST

# Blood Pressure Monitor

- 5-Year Limited Warranty

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---

## **INTRODUCTORY INFORMATION:**

<b>Important Product Notices and Safety Instructions .....</b>	<b>3</b>
<b>About Blood Pressure .....</b>	<b>5</b>
<b>Blood Pressure Standard.....</b>	<b>6</b>
<b>How This Blood Pressure Monitor Works.....</b>	<b>8</b>

## **USAGE INFORMATION:**

<b>Name/Function of Each Part .....</b>	<b>9</b>
<b>Display Explanations .....</b>	<b>10</b>
<b>Installing Batteries .....</b>	<b>12</b>
<b>Wrist Positioning Guide .....</b>	<b>14</b>
<b>Wrist Positioning Guide Settings .....</b>	<b>15</b>
<b>Date &amp; Time Set Procedure.....</b>	<b>16</b>
<b>Bluetooth® Operation .....</b>	<b>19</b>
<b>Applying the Cuff.....</b>	<b>20</b>
<b>Correct Measurement Posture .....</b>	<b>22</b>

<b>Measurement Procedure.....</b>	<b>23</b>
<b>Transmit Readings to Your Mobile Device... </b>	<b>27</b>
<b>Earn Balance® Reward Points .....</b>	<b>29</b>
<b>Risk Category Index .....</b>	<b>30</b>
<b>Irregular Heartbeat Detector (IHB) .....</b>	<b>31</b>
<b>Recalling Values from Memory.....</b>	<b>32</b>
<b>Clearing Values from Memory .....</b>	<b>33</b>
<b>Important Notes Regarding Your Blood Pressure Measurement.....</b>	<b>34</b>
<b>Care, Maintenance &amp; Cleaning .....</b>	<b>35</b>
<b>Potential for Electromagnetic Interference.....</b>	<b>37</b>
<b>Electromagnetic Compatibility (EMC).....</b>	<b>38</b>
<b>Troubleshooting .....</b>	<b>40</b>
<b>Specifications .....</b>	<b>42</b>
<b>5 Year Limited Warranty .....</b>	<b>44</b>

# IMPORTANT PRODUCT NOTICES AND SAFETY INSTRUCTIONS

---

When using your blood pressure monitor, basic precautions should always be followed. Please read and follow all instructions and warnings before using this product. Save these instructions for future reference.

- Please note that this is a home healthcare product only and it is not intended to serve as a substitute for the advice of a physician or medical professional.
- This device uses the oscillometric method to measure systolic and diastolic blood pressure, as well as heart rate.
- Do not use this device for diagnosis or treatment of any health problem or disease. Measurement results are for reference only. Consult a healthcare professional for interpretation of blood pressure measurements. Contact your physician if you have or suspect any medical problem. Do not change your medications without the advice of your physician or healthcare professional.
- This product is not suitable for people with arrhythmias.
- This device may have difficulty determining the proper blood pressure for pregnant women and for individuals with irregular heartbeat, diabetes, poor circulation of blood, kidney problems or for users who have suffered from a stroke.
- Excessive use may result in blood flow interference, which is likely to cause uncomfortable sensations, such as partial subcutaneous hemorrhage, or temporary numbness to your wrist. In general these symptoms should not last long. However, if you do not recover in time, please seek advice from a medical professional.

- Consult your physician before measuring blood pressure at the wrist if you have any of the following conditions: Severe hypertension, severe diabetes, severe arteriosclerosis, severe kidney disease or vascular disease that may compromise circulation.
- The pulse display is not suitable for checking the frequency of heart pacemakers.
- Electromagnetic interference: The device contains sensitive electronic components. Avoid strong electrical or electromagnetic fields in the direct vicinity of the device (e.g. mobile telephones, microwave ovens). These may lead to temporary impairment of measurement accuracy.
- Consider the electromagnetic compatibility of the device (ex. power disturbance, radio frequency interference etc.) Please use this device indoors only.
- Use blood pressure monitor only for its intended use.
- Do not wrap the cuff around body parts other than your wrist.
- Not for use by or on persons under the age of 18.
- Do not use this device on infants, children, or those who cannot express their own intention.
- Use only 1.5V “AAA” alkaline batteries for power supply.

Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultatory method, within the accuracy limits prescribed by the American National Standard, *Manual, electronic, or automated sphygmomanometers*.

# ABOUT BLOOD PRESSURE

---

## ***What is Blood Pressure?***

Blood pressure is the pressure exerted on the artery walls while blood flows through the arteries. The pressure measured when the heart contracts and sends blood out of the heart is systolic (highest) blood pressure. The pressure measured when the heart dilates with blood flowing back into the heart is called diastolic (lowest) blood pressure.

## ***Why Measure Your Blood Pressure?***

Among today's various health problems, those associated with high blood pressure are very common. High blood pressure dangerously correlates with cardiovascular diseases. Therefore, blood pressure monitoring is important for identifying those at risk.

## ***Why Do My Readings Vary?***

Blood pressure is a body parameter that is subject to normal variations throughout the day. A single reading that is different from yours or your doctor's readings are not necessarily inaccurate. The average of several readings, taken under similar conditions, using the same arm is preferred for accurate blood pressure readings.

## ***Why Are My Readings Different Than Those Taken at My Doctor's Office?***

Many experience a phenomenon called "White Coat Hypertension" when measured by a doctor. White Coat Hypertension refers to blood pressure that rises above its usual level when measured in a clinical setting, such as a doctor's office.

# BLOOD PRESSURE STANDARD

---

The table below contains defined levels for hypertension that are publicly available from the National Heart Lung and Blood Institute at the U.S. National Institutes of Health (NIH) ([http://www.nhlbi.nih.gov/health/dci/Diseases/Hbp/HBP\\_WhatIs.html](http://www.nhlbi.nih.gov/health/dci/Diseases/Hbp/HBP_WhatIs.html)). Users can compare their own blood pressure readings against these defined levels to determine if they may be potentially at increased risk.

This table is applicable to most adults aged 18 and older.

Category	Systolic (mmHg)		Diastolic (mmHg)
Normal	Less than 120	and	Less than 80
Prehypertension	120-139	or	80-89
High Blood Pressure			
Stage 1	140-159	or	90-99
Stage 2	160 or higher	or	100 or higher

Blood pressure tends to go up and down, even in people who normally don't have high readings. If your numbers stay above the "normal" range most of the time, you may be at increased risk and should consult your physician.

Although one can easily find where their own blood pressure readings fall on this table, this monitor comes equipped with a Risk Category Index that automatically compares each reading to the defined levels and provides a helpful cue if your reading falls into one of the stages that could potentially indicate increased risk. See page 30 for more information on this feature.

Please note that cues provided by this monitor are only intended to assist you in using this table. The table and cues are only provided for convenience to help you understand your non-invasive blood pressure reading as it relates to the NIH information. They are not a substitute for a medical examination by your physician. It is important for you to consult with your physician regularly. Your physician will tell you your normal blood pressure range as well as the point at which you may actually be considered to be at risk.

## HOW THIS BLOOD PRESSURE MONITOR WORKS

---

This monitor uses innovative GentleInflate™ technology to detect your blood pressure. With one touch of a button, the monitor automatically takes the reading during inflation, then deflates, delivering a quick, comfortable and quiet measurement every time.

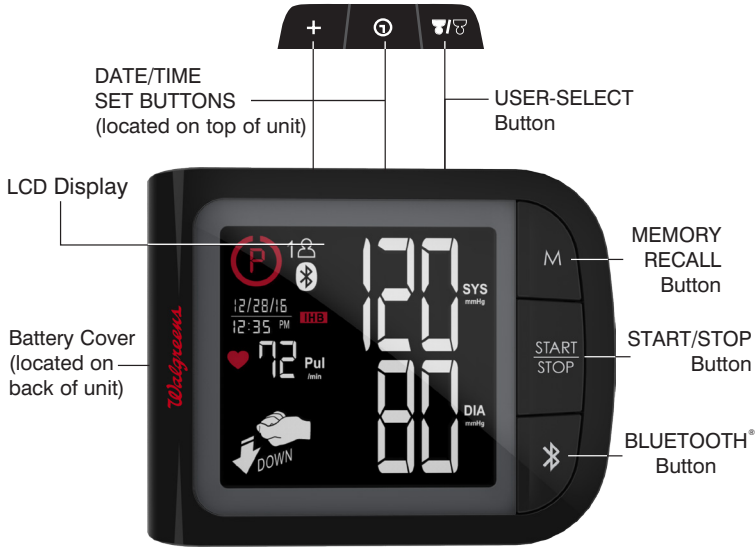
To help guide your wrist into the ideal position, this blood pressure monitor comes equipped with an optional wrist positioning guide to help you take an accurate blood pressure reading. Please see page 14 for more information on this feature.

Please note that any muscle movement during inflation will cause measurement error. When measurement is complete, the monitor will display your systolic pressure, diastolic pressure, and pulse readings.


The monitor automatically finds where your measurement results fall on the NIH's National Heart Lung and Blood Institute's table and provides a cue if your reading falls into one of the stages that could potentially indicate increased risk. See page 30 for more information on this feature.

The appearance of the **IHB** icon indicates that a pulse irregularity consistent with an irregular heartbeat was detected during measurement. Refer to page 31 for more information on the Irregular Heartbeat Detector.

# NAME/FUNCTION OF EACH PART


## Accessories for unit :

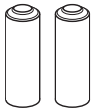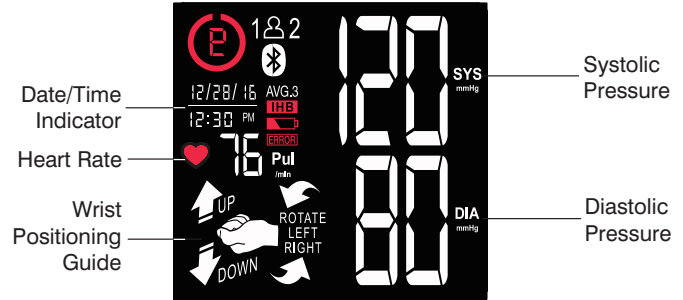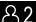
2 "AAA" 1.5V alkaline batteries and a carrying case are included with the unit.

# DISPLAY EXPLANATIONS

## Display Symbols:

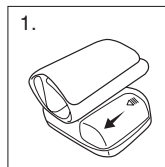
	<b>User 1:</b> Appears when the monitor is operated by User 1.
	<b>User 2:</b> Appears when the monitor is operated by User 2.
	<b>Weak Battery Symbol:</b> Appears when batteries should be replaced.
	<b>Pulse Symbol:</b> Shows the heart rate per minute.
	<b>Risk Category Index:</b> See page 30 for more information.
	<b>Irregular Heartbeat Detector:</b> See page 31 for more information.
<b>AVG. 3</b>	<b>Memory Average:</b> Displays average of last 3 readings.
	<b>Wrist Positioning Guide:</b> See page 14 for more information.
	<b>Data Transmitting Symbol:</b> Appears when data is transmitting to your mobile device.

If **ERROR** and any of the following letters and numbers appear in the area that systolic pressure should be displayed or the wrist icon appears, an error has occurred with your reading. See Troubleshooting section of this manual for more information.


<b>EE</b>	<b>Measurement Error:</b> Wrap the cuff correctly and keep wrist steady during measurement. Measure again.
<b>E1</b>	<b>Air Circuit Abnormality:</b> Measure again.
<b>E2</b>	<b>Pressure Exceeding 300 mmHg:</b> Switch the unit off to clear, then measure again.
<b>E3</b>	<b>Error Determine Blood Pressure Measurement:</b> Rewrap cuff properly, keep steady, and measure again.
<b>E4</b>	<b>Data Transmission Error:</b> Monitor cannot connect to the mobile device to transmit data. Make sure <i>Bluetooth</i> is ON.
<b>EP</b>	<b>System Error:</b> Measure again.

## INSTALLING BATTERIES

1. The battery cover is located on the back of the monitor. Remove battery cover by pressing down and pulling away from the monitor.

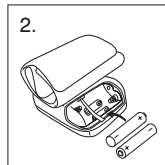
2. Install batteries according to the polarity indications inside the compartment (as shown in Figure 2). **Battery Type: 2 Alkaline LR03 (AAA) size.**


**Note:** For easy battery removal, place batteries on top of the ribbon located inside of the battery compartment.

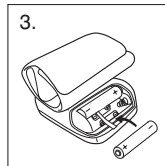3. Insert the other battery, according to the polarity indications inside the compartment (as shown in Figure 3).

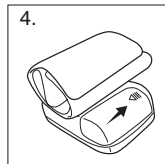4. Replace the cover by inserting it back into place.

### Replace the batteries if:

- The weak battery symbol appears.
- When any button is pressed and nothing is displayed on the screen.

### Note:

- If batteries are removed or replaced, the date and time will need to be re-set either manually or automatically using your *Bluetooth* mobile device.
- Replace all batteries at one time (as simultaneous set). Use only 1.5V “AAA” alkaline batteries.
- When the batteries are removed, the measurement values stored in memory are retained. However, the date and time must be reset.
- Remove batteries when unit is not in use for extended periods of time.
- Clean contacts on battery and in battery compartment with a soft dry cloth each time you install batteries.
- **Batteries are hazardous waste. Do not dispose of them together with the household garbage.**

## WRIST POSITIONING GUIDE

- When using wrist blood pressure monitors, it is important to take the measurement while relaxing your wrist at heart level. This will help ensure an accurate reading will be taken. Because it is difficult to find the ideal location for a wrist measurement, we have included an optional Wrist Positioning Sensor that will assist in directing your wrist to the ideal location for a blood pressure measurement.
- After you press **START/STOP button**, the display will illuminate with different icons that are designed to help you move your wrist. Once the ideal location is found, the Pulse Symbol will flash and measurement will begin. Keep your arm still until measurement is completed.

SYMBOL	ACTION
	Move your wrist up.
	Move your wrist down.
	Rotate your wrist to your left.
	Rotate your wrist to your right.

### Note:

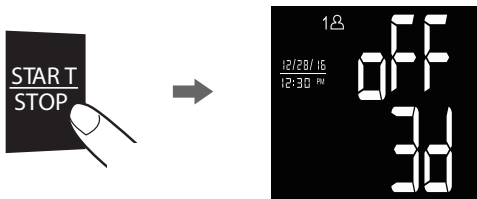
- Even if the monitor is not in the ideal position after 10 seconds, the measurement will still begin. If this occurs, the wrist symbol  will appear and be stored in the memory with the measurement.


# WRIST POSITIONING GUIDE SETTINGS

**Note:** By default, the Wrist Positioning Guide is ON.

1. To turn the Wrist Positioning Guide OFF, press and hold the **START/STOP button** for 3 seconds.

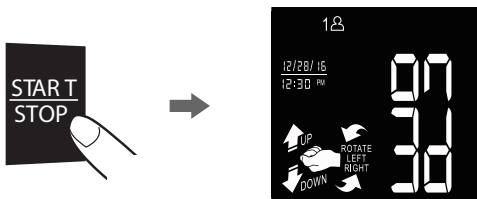2. To turn the Wrist Positioning Guide ON, press and hold the **START/STOP button** again for 3 seconds.

# DATE & TIME SET PROCEDURE

Date and time can be set by two methods; either sync automatically using your *Bluetooth* mobile device, or set manually using the **SET buttons (⊖ & +)** on the monitor.

If setting the date and time by syncing with your mobile device, it is important that this is done prior to taking any measurements to ensure the date and time are accurate.

**To set using your Bluetooth® mobile device:**

1. Before taking the first measurement, open the **Walgreens Connect App** on your device. If you have not already installed the **Walgreens Connect App** onto your device, it is available to download on the App Store<sup>SM</sup> and on Google Play<sup>TM</sup>. Make sure you have the app downloaded and open on your mobile device before trying to sync the date and time.

2. Press the **Bluetooth button**  on the monitor to enable the connection.

3. Once the *Bluetooth* connection is established, the date and time will automatically update on the blood pressure monitor.

### Note:

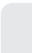
- Once the date and time have been successfully synced, future readings will automatically have the correct date and time.
- If you are having trouble automatically syncing the date and time, make sure *Bluetooth* is ON on your device and ON on the blood pressure monitor.


### To set manually:

1. To set the date and time, press the **DATE/TIME SET button** .

2. The display will show a blinking number showing the MONTH. Change the MONTH by pressing the **+ button**. Each press will increase the number by one in a cycling manner. Press the **DATE/TIME SET button**  again to confirm the entry, and the screen will show a blinking number representing the DAY.

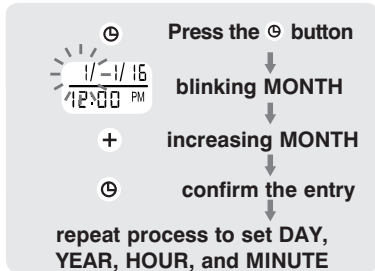3. Change the DAY, YEAR, HOUR, & MINUTE as described in Step 2 above, using the **+ button** to change the numbers and the **DATE/TIME SET button**  to confirm the entries.

### Note:

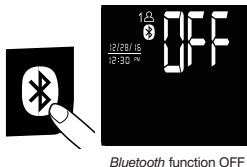
- The date and time will only need to be set manually if the monitor will not be used with the **Walgreens Connect App**.
- If the **Walgreens Connect App** is used after manually setting the date and time on the blood pressure monitor, the date and time on the mobile device will override the date and time on the blood pressure monitor.

# BLUETOOTH® OPERATION

This monitor has *Bluetooth* function ON for your convenience. This will allow your readings to automatically transmit to the **Walgreens Connect App**. By turning this function OFF, the measurements cannot be transmitted.

## Turning Bluetooth® Function OFF:


While the screen is off (nothing is shown on the screen), press and hold the **Bluetooth button**  for 3 seconds to turn *Bluetooth* OFF.

- When *Bluetooth* is OFF, measurements cannot be transmitted to the **Walgreens Connect App**.

## Turning Bluetooth® Function ON:

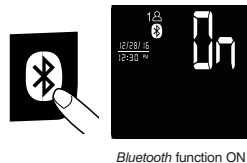While the screen is off (nothing is shown on the screen), press and hold the **Bluetooth button**  for 3 seconds to turn *Bluetooth* ON.


- When *Bluetooth* is ON, readings can automatically and manually be transmitted to the App.
- The **Walgreens Connect App** must be open on your mobile device in order to establish a connection with the monitor.

Bluetooth QDID Profile Information:  
QDID: 45008, 51582, & 52727


# APPLYING THE CUFF

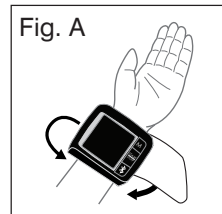1. Remove all watches, wrist jewelry, etc. prior to attaching the wrist monitor. Clothing sleeves should be rolled up and the cuff should be wrapped around bare skin for correct measurements.


2. Apply cuff to left wrist with palm facing up as shown in Fig. A.

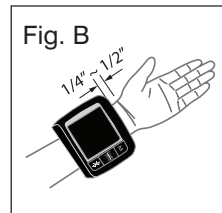
3. Make sure the edge of the cuff is about 1/4"~1/2" (1 cm) from the palm as shown in Fig. B.

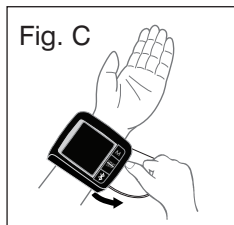4. In order to ensure accurate measurements, fasten the hook and loop strap securely around your wrist so there is no extra space between the cuff and the wrist as shown in Fig. C. If the cuff is not wrapped tight enough, the measurement values will not be accurate.

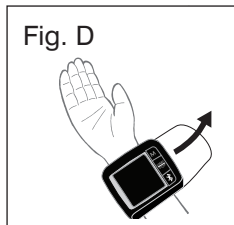5. If your physician has diagnosed you with poor circulation in your left arm, place the cuff around your right wrist as shown in Fig. D.

**Note:**

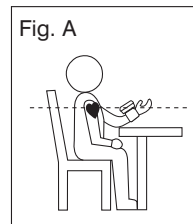
- This device should not be used when your wrist is wounded or injured.
- If it is not possible to use the cuff on your left wrist, it can also be used on your right wrist. However, all measurements should be made using the same wrist.

## CORRECT MEASUREMENT POSTURE


1. Place your elbow on a table so that the cuff is at the same level as your heart as shown in Fig. A. Relax your entire body, especially between your elbow and fingers.

**Note:**

Your heart is located slightly below your armpit.

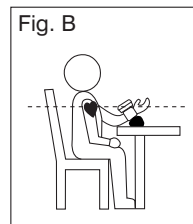
2. If the cuff is not at the same level as your heart or if you cannot keep your arm completely still throughout the reading, use a soft object such as a folded towel to support your wrist as shown in Fig. B.

3. Turn your palm upwards.

4. Sit upright in a chair, and take 5-6 deep breaths.

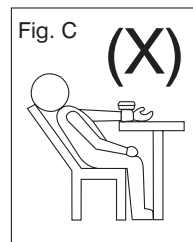Avoid leaning back while the measurement is being taken as shown in Fig. C.

# MEASUREMENT PROCEDURE

## Note:

- Position the monitor at the same level as your heart during measurement to ensure accurate readings.
- Blood pressure changes with every heartbeat and is in constant fluctuation throughout the day.
- Blood pressure measurement can be affected by the position of the user, his or her physiologic condition and other factors.

## Before Measurement:

To help ensure an accurate reading, follow these instructions before taking a measurement:

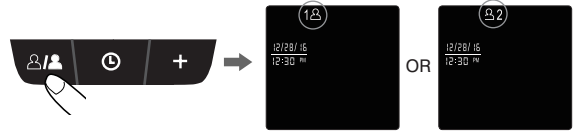
- Wait 1 hour after exercising, bathing, eating, drinking beverages with alcohol or caffeine, or smoking, to measure blood pressure.
- Sit quietly and rest for 15 minutes.
- Stress raises blood pressure. Avoid taking measurements during stressful times.
- Take your reading in a comfortable environment as measurements can be affected by hot or cold temperatures. Take your blood pressure at normal body temperature.

## During Measurement:

- Do not talk or move your arm or hand muscles.
- Do not cross your legs. Sit with feet flat on the floor.
- Dot not touch cuff or monitor during measurement.

If you are using this blood pressure monitor for the first time, please remove the protective film from the screen.

1. Press the **USER-SELECT** button to choose User 1 or User 2.


2. With the cuff wrapped around your wrist, press the **START/STOP** button. Do not start the measurement process unless the cuff is wrapped around your wrist.
3. If your wrist is not in the ideal position, the screen will suggest moving your wrist up or down, or rotating to the left, or to the right.

SYMBOL	ACTION
	Move your wrist up.
	Move your wrist down.
	Rotate your wrist to your left.
	Rotate your wrist to your right.


**Note:**

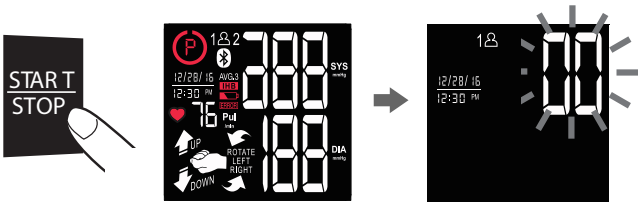When the Wrist Positioning Guide Function is turned OFF, this instruction will not appear.

- All digits will light up, checking the display functions. The checking procedure will be completed after about 3 seconds.

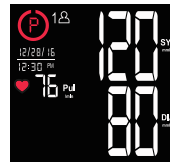- After all symbols disappear, the display will show “00”. The monitor is “**Ready to Measure**” and will **automatically** inflate the cuff to begin measurement.
- When the measurement is completed, the cuff will deflate entirely and systolic pressure, diastolic pressure and pulse will be shown simultaneously on the LCD screen. **The measurement is then automatically stored into memory.**

**Note:**

- This monitor automatically turns off approximately 1 **minute** after last operation. You may also press the **START/STOP button** to turn the unit off.
- To interrupt the measurement, you may press the **START/STOP button**. The cuff will deflate immediately after this button is pressed.

# TRANSMIT READINGS TO YOUR MOBILE DEVICE

This monitor syncs your blood pressure readings to the **Walgreens Connect App** that is free to download on the App Store<sup>SM</sup> and on Google Play<sup>TM</sup>. Make sure you have the **Walgreens Connect App** downloaded, and open on your mobile device before trying to transmit your blood pressure measurements.

Before attempting to sync the blood pressure monitor with your mobile device, make sure *Bluetooth* is turned ON, on both your mobile device and the monitor.

## Automatically transmit readings:

After a measurement is taken, the *Bluetooth* icon  will appear on the screen as the monitor automatically transmits your blood pressure readings to the App.

## Manually transmit readings:

Press the **Bluetooth button** . The readings will automatically be transmitted to the App. The *Bluetooth* icon  will appear on the screen as the monitor transmits your blood pressure readings to the App.

If the transmission(s) were successful, the *Bluetooth* icon  will be displayed on the screen. See Fig. 1.

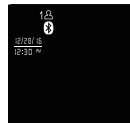
Fig. 1

If the transmission(s) were unsuccessful, E4 and the **ERROR** symbol will be displayed.

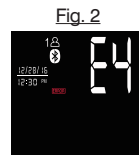
Fig. 2

## Note:

- The measurement is stored in the monitor's memory even if the reading is not transmitted to your mobile device.
- **The Walgreens Connect App** must be open on your mobile device in order to transmit your measurements.
- Only new readings will be accepted by the App.
- The blood pressure monitor can only be assigned to one Walgreens user account.
- To ensure readings transmitted to the App have the correct date and time, it is important that the correct date and time is set on the blood pressure monitor before taking measurements.
- Measurements transmitted to the App cannot be edited.
- For a complete list of compatible mobile devices please visit [www.walgreens.com/connectapp](http://www.walgreens.com/connectapp).

## EARN BALANCE® REWARDS POINTS\*

Earn Balance® Rewards points for taking blood pressure readings when you sync this device to your smartphone and use the **Walgreens Connect App**.

\*Balance® Rewards:

- Limitations apply to earning points.
- Limit 20 points per blood pressure test; one log per day.
- Information provided to Walgreens online is covered by the terms of our Online Privacy and Security Policy and the terms and conditions of Balance® Rewards.
- Personally identifiable information reported to earn points is not covered under HIPAA or the Walgreens Notice of Privacy Practices.
- See full terms and details at [www.walgreens.com/healthychoices](http://www.walgreens.com/healthychoices).

## RISK CATEGORY INDEX

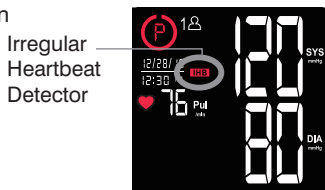
This monitor comes equipped with a Risk Category Index that automatically compares each reading to defined levels established by the U.S. National Institutes of Health's (NIH) National Heart Lung and Blood Institute as described earlier in this manual, and provides a helpful cue if your reading falls into one of the stages that could potentially indicate increased risk. No cue is given if the reading falls in the normal range as defined by NIH. Please note that cues provided by this monitor are only intended to assist you in using this table. The table and cues are only provided for convenience to help you understand your non-invasive blood pressure reading as it relates to the NIH information. They are not a substitute for a medical examination by your physician. It is important for you to consult with your physician regularly. Your physician will tell you your normal blood pressure range as well as the point at which you may actually be considered to be at risk.

	Stage 2 Hypertension		Prehypertension
	Stage 1 Hypertension		Normal (No Symbol)

# IRREGULAR HEARTBEAT DETECTOR (IHB)

The appearance of the **IHB** icon indicates that a pulse irregularity consistent with an irregular heartbeat was detected during measurement. Usually this is not a cause for concern. However, if the symbol appears often, we recommend you seek medical advice. Please note that the device does not replace a cardiac examination, but serves to detect pulse irregularities at an early stage.

Irregular Heartbeat Detector

Movement, shaking or talking during the measurement can result in pulse irregularities that may cause the appearance of this icon. Therefore, it is of great importance to not move or talk during measurement.

To determine the presence of an irregular heartbeat, the average of the heartbeat intervals is calculated with the first 3 normal effective heartbeat values. It is important to note that the average is not a strict mathematical averaging of all recorded intervals. At least 3 beats with 25% or greater difference from the average heartbeat interval will generate the **IHB** icon on the screen.

## IMPORTANT INFORMATION:

This blood pressure monitor is not designed for use by people with arrhythmias nor for diagnosing or treating an arrhythmia problem. As a safeguard, we recommend that if you have arrhythmias such as atrial or ventricular premature beats and atrial fibrillation or any other special conditions you should check with your physician before using your blood pressure monitor.

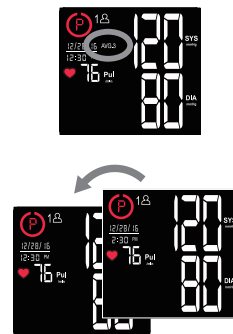
# RECALLING VALUES FROM MEMORY

This monitor can be used by two individuals. **Each user can store up to 120 measurements.**

1. Press the **USER-SELECT** button on top of the monitor and choose User 1 or User 2.

2. Press the **M** button to access the memory.
3. The monitor will first display the calculated average applied to the last 3 memories ("AVG. 3").
4. Every new press of the **M** button will recall a previous reading. The latest reading will be recalled first.


## Note:


If the wrist positioning guide is ON before the measurement and the monitor is not in the ideal position, the  icon will be displayed with the saved readings.

## CLEARING VALUES FROM MEMORY

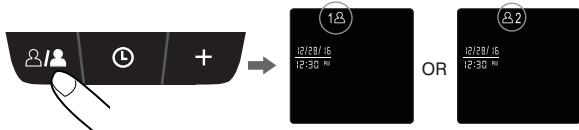1. Press the **USER-SELECT** button on top of the monitor to select User 1 or User 2.

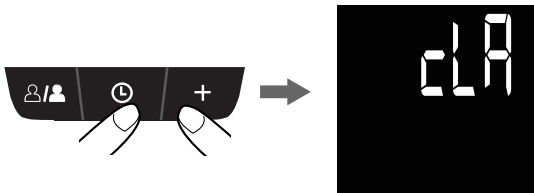2. Press **M** button to enter memory recall mode.
3. Press and hold the **DATE/TIME SET**  $\ominus$  and **+** buttons at the same time and the data for the selected user will be erased automatically.

### Note:

Once deleted, your readings cannot be restored.

## IMPORTANT NOTES REGARDING YOUR BLOOD PRESSURE MEASUREMENT

- It is suggested that you take your measurements at the same time each day and use the same wrist for consistency.
- Users should wait a minimum of 5 minutes before taking additional measurements. More time may be necessary depending upon your physiology.
- The measurement results that users receive are for reference only. If users have any blood pressure concerns, please consult a physician.
- Once inflation reaches 300 mmHg, the unit will deflate automatically for safety reasons.
- This product is not suitable for people with arrhythmias.
- This device may have difficulty determining the proper blood pressure for pregnant women and for users with irregular heartbeat, diabetes, poor circulation of blood, kidney problems, or for users who have suffered a stroke.

# CARE, MAINTENANCE & CLEANING

---

- Clean the blood pressure monitor body and cuff carefully with a slightly damp, soft cloth. Do not press. Do not wash cuff or use chemical cleaner on it. Never use thinner, alcohol or petrol (gasoline) as cleaner.
- Leaky batteries can damage the unit. Remove the batteries when the unit will not be used for a long time.
- Follow local ordinances and recycling instructions regarding disposal or recycling of the device and device components, including batteries.
- If the unit is stored near freezing, allow it to acclimate to room temperature before use.
- This blood pressure monitor is not field serviceable. You should not use any tool to open the device nor should you attempt to adjust anything inside the device. If you have any problems with this device, please contact HoMedics Consumer Relations (contact information on warranty page).
- Do not immerse the unit in water as this will result in damage to the unit.
- Do not subject the monitor or cuff to extreme temperatures, humidity, moisture, or direct sunlight. Protect from dust.

- Do not fold the cuff tightly.
- Do not disassemble the monitor or cuff. If in need of repair, refer to the warranty section of this manual.
- Do not subject the monitor to extreme shocks (do not drop on floor).
- Do not inflate the cuff unless wrapped around wrist.
- Do not wrap the cuff around body parts other than your wrist.
- Do not drop or insert any object into any opening.
- This monitor may not meet its performance specifications if stored or used outside of these temperature and humidity ranges:

- **Storage/Transportation Environment**

Temperature:  
-13°F ~ 158°F  
(-25°C ~ 70°C)

Humidity:  
Less than 93% RH

- **Operation Environment**

Temperature:  
41°F ~ 104°F  
(5°C ~ 40°C)

Humidity:  
15 ~ 93% RH

# POTENTIAL FOR ELECTROMAGNETIC INTERFERENCE

To avoid inaccurate results caused by electromagnetic interference between electrical and electronic equipment, do not use the device near a cell phone or microwave oven. For most wireless communication devices, it is recommended to maintain a distance of 10.8 feet (3.3m) in order to avoid electromagnetic interference.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**FEDERAL COMMUNICATIONS COMMISSION COMPLIANCE STATEMENT**  
Changes or modifications to this equipment not expressly approved by the manufacturer could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

# ELECTROMAGNETIC COMPATIBILITY (EMC)

- **Guidance and manufacturer's declaration - electromagnetic emissions**  
The device is intended for use in the electromagnetic environments listed below, and should only be used in such environments:

Emissions test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	RF energy is used only to maintain device's operation. Therefore, its RF emissions are so low that it's not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Not Applicable	The device is suitable for use in all establishments, including domestic establishments, and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Voltage fluctuations/flicker emissions IEC 61000-3-3	Not Applicable	

- **Guidance and manufacturer's declaration - electromagnetic immunity**  
The device is intended for use in the electromagnetic environments listed below, and should only be used in such environments:

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 8 kV contact ± 15 kV air	The relative humidity should be at least 5 %.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment. Recommended separation distance $r = \frac{l}{188}$ (m) where l is the current in amperes in a power bus or an appliance wire and r is the recommended separation distance between your device and the power bus or appliance wire, in meters (m).

- **Recommended separation distances between portable and mobile RF communication equipment and the device.**  
The device is intended for use in an electromagnetic environment where radiated RF disturbances are under control. User can help prevent electromagnetic interference by keeping the device at a minimum distance from portable and mobile RF communications equipment (transmitters). Below table details the maximum output power of transmitter:

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz Not Applicable	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$
0.01	NA	0.12	0.23
0.1	NA	0.38	0.73
1	NA	1.2	2.3
10	NA	3.8	7.3
100	NA	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

• **Guidance and manufacturer's declaration – electromagnetic immunity**

The device is intended for use in the electromagnetic environments listed below, and should only be used in such environments:

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	Not Applicable	Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.3 \sqrt{P}$ 800 MHz to 2.5 GHz  where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).  Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, a should be less than the compliance level in each frequency range. b Interference may occur in the vicinity of equipment marked with the following symbol: 
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	10 V/m 80 MHz to 2.7 GHz	

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the device.

b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.

# TROUBLESHOOTING

If any abnormality arises during use, please check the following points:

SYMPTOMS	POSSIBLE CAUSES	CORRECTION
Unit does not turn on when the <b>START/STOP button</b> is pushed.	Batteries have run down.	Replace them with two new "AAA" alkaline batteries.
	Battery polarities have been positioned incorrectly.	Re-insert the batteries in the correct positions.
EE measurement error symbol shown on display or the blood pressure value is displayed excessively low (or high).	Cuff has been placed incorrectly.	Wrap the cuff properly so that it is positioned correctly. Measure again.
	Did you talk or move during measurement?	Keep wrist steady during measurement. Measure again.
	Shaking of the wrist with the cuff on.	
E1 error symbol shown on display	Air circuit abnormality.	Measure again.
E2 error symbol shown on display	Inflation pressure exceeding 300 mmHg.	Switch the unit OFF, then measure again.
E3 error symbol shown on display	Can't determine blood pressure measurement data.	Rewrap the cuff properly, keep steady, and measure again.

**Note:** If the unit still does not work, contact HoMedics Consumer Relations. Under no circumstance should you disassemble or attempt to repair the unit by yourself. Contact information for HoMedics Consumer Relations Department can be found on the warranty page.

E4 - Data transmission error symbol shown on the display	<i>Bluetooth</i> function is turned OFF, on your mobile device.	Turn <i>Bluetooth</i> ON, on your mobile device, and try again.
	<i>Bluetooth</i> function is turned OFF, on the blood pressure monitor.	Turn the <i>Bluetooth</i> function ON, on your blood pressure monitor (see page 19), then try again.
	Mobile device does not pair with the blood pressure monitor.	Make sure <i>Bluetooth</i> is ON, on the mobile device and the blood pressure monitor, and try again.
		Requires compatible mobile device. This blood pressure monitor is compatible with mobile devices that are running iOS 5.0 or newer, or Android™ 4.3 and newer. For a complete list, go to <a href="http://www.walgreens.com/connectapp">www.walgreens.com/connectapp</a> .
	The App on the mobile device is not open.	Make sure the app is open on the mobile device, and try again.
	The blood pressure monitor and mobile device are out of transmitting range.	Make sure the mobile device and blood pressure monitor are within the acceptable range of 32 feet (10 meters).
	Unexpected loss of electrical/mechanical integrity.	Make sure your last reading is stored in memory and the App is open and try again. See page 32.
		Remove batteries, re-insert, and try again.
Return the device to your local distributor or importer.		
EP error symbol shown on display	System error.	Measure again.

## SPECIFICATIONS

<b>Power Source:</b>	(2)1.5V DC “AAA” Alkaline batteries
<b>Measurement Method:</b>	Oscillometric
<b>Measurement Range:</b>	Pressure: 40~280 millimeters Mercury (mmHg) Pulse: 40~199 beats/minute
<b>Accuracy:</b>	Pressure: ±3 mmHg Pulse: ±5% of reading
<b>Pressure Sensor:</b>	Semi-conductor
<b>Inflation:</b>	Pump driven
<b>Deflation:</b>	Automatic Pressure release valve
<b>Memory Capacity:</b>	120 memories for each User (240 total)
<b>Auto-shut-off:</b>	1 minute after last button operation
<b>Operation Environment:</b>	Temperature: 41°F~104°F (5°C~40°C) Humidity: 15% ~ 93% RH
<b>Storage/Transportation Environment:</b>	Temperature: -13°F~158°F (-25°C~70°C) Humidity: Less than 93% RH

<b>Weight:</b>	0.21 lbs (96 g) (without batteries)
<b>Wrist Circumference:</b>	5.3" - 7.7" (13.5cm - 19.5cm)
<b>Accessories:</b>	(2) "AAA" alkaline batteries, instruction manual, hard plastic storage case
<b>Dimensions:</b>	3.58" (L) x 1.42" (W) x 2.76" (H) 91mm(L) x 36mm(W) x 70mm(H)
	Internally Powered equipment Type BF Equipment IP22-Degrees of protection provided by enclosures Not suitable for use in presence of flammable anesthetic mixture with air, oxygen or nitrous oxide. Continuous operation with short-time loading.
	Follow instructions for use.

**Note:** These specifications are subject to change without notice.

Made exclusively for

*Walgreens*

**For service or repair, do not return this unit to the retailer. Contact HoMedics Consumer Relations at:**

**Email:**  
**cservice@**  
**homedics.com**

**Phone:**  
**1-800-466-3342**  
**Business Hours:**  
**8:30am-7pm ET**  
**Monday-Friday**

**Distributed by**  
**HOMEDICS**<sup>®</sup>  
HoMedics, LLC.  
3000 Pontiac Trail  
Commerce  
Township, MI  
48390

Printed in China

Learn more at:  
[www.walgreensbp.com](http://www.walgreensbp.com)

LIMITED FIVE YEAR WARRANTY

HoMedics sells its products with the intent that they are free of defects in manufacture and workmanship for a period of five years from the date of original purchase, except as noted below. HoMedics warrants that its product will be free of defects in material and workmanship under normal use and service. This blood pressure monitor meets the simulated measurement cycles test requirement per EN1060-3, part 8.10. This warranty extends only to consumers and does not extend to Retailers.

To obtain warranty service on your HoMedics product, contact a Consumer Relations Representative by telephone at 1-800-466-3342 for assistance. Please make sure to have the model number of the product available.

HoMedics does not authorize anyone, including, but not limited to, Retailers, the subsequent consumer purchaser of the product from a Retailer or remote purchasers, to obligate HoMedics in any way beyond the terms set forth herein. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; improper installation; unauthorized repairs or modifications; improper use of electrical/power supply; loss of power; dropped product; malfunction or damage of an operating part from failure to provide manufacturers recommended maintenance; transportation damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of HoMedics.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adoption to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. THERE SHALL BE NO OTHER WARRANTIES EXPRESS OR IMPLIED INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS OR ANY OTHER OBLIGATION ON THE PART OF THE COMPANY WITH RESPECT TO PRODUCTS COVERED BY THIS WARRANTY. HOMEDICS SHALL HAVE NO LIABILITY FOR ANY INCIDENTAL, CONSEQUENTIAL OR SPECIAL DAMAGES. IN NO EVENT SHALL THIS WARRANTY REQUIRE MORE THAN THE REPAIR OR REPLACEMENT OF ANY PART OR PARTS WHICH ARE FOUND TO BE DEFECTIVE WITHIN THE EFFECTIVE PERIOD OF THE WARRANTY. NO REFUNDS WILL BE GIVEN. IF REPLACEMENT PARTS FOR DEFECTIVE MATERIALS ARE NOT AVAILABLE, HOMEDICS RESERVES THE RIGHT TO MAKE PRODUCT SUBSTITUTIONS IN LIEU OF REPAIR OR REPLACEMENT.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or sales of such products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate as to any products or parts thereof which are repaired, replaced, altered, or modified, without the prior express and written consent of HoMedics.

This warranty provides you with specific legal rights. You may have additional rights which may vary from state to state. Because of individual state regulations, some of the above limitations and exclusions may not apply to you.

**The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Walgreen Company is under license. Other trademarks and trade names are those of their respective owners.**

**Android™ and Google Play™ are trademarks of Google, Inc.**

**HoMedics is a registered trademark of HoMedics, LLC.  
GentleInflate™ is a trademark of HoMedics, LLC.  
© 2016 HoMedics, LLC. All rights reserved.**

**App Store™ is a service mark of Apple Inc.**

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**Android™ y Google Play™ son marcas comerciales de Google, Inc.**
**HoMedics® es una marca registrada de HoMedics, LLC.
GentleInflate™ es una marca registrada de HoMedics, LLC.
© 2016 HoMedics, LLC. Todos los derechos reservados.**
**App Store™ es una marca de servicio de Apple Inc.**